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Tab. 7

Table of Contents

1 Explanation of Symbols	39	7 Handling the Coolant/Lubricant (C/L)	62
2 Safety Instructions	40	7.1 C/L analysis	62
2.1 General Safety Instructions	40	7.2 Removing the tank from inside the cart	63
2.2 Unattended operation	41	7.3 Emptying the C/L container	63
2.3 Processing titanium blanks	41	7.4 Cleaning the tank	64
3 Suitable Personnel	42	7.5 Filling in coolant/lubricant	64
4 Product Specifications	42	7.6 Installing the tank into the cart	64
4.1 Processing machine	42	7.7 Cleaning the C/L system	65
4.2 Accessories	48	7.8 Changing the filter cartridge	65
5 Installation	51	8 Cleaning and Maintenance	65
5.1 Processing machine	51	8.1 Processing machine	65
5.2 Setting up the tank	53	8.2 Tank	68
5.3 Service unit	53	9 Malfunctions, Repairs and Warranty	69
5.4 Installation of the CAM and machine software	55	9.1 Malfunctions	69
5.5 Connecting the interfaces	56	9.2 Repairs	69
5.6 Connection to the PC	56	9.3 Warranty	69
5.7 Establishing the network connection	57	10 Environmental Protection	69
5.8 Calibration of the processing machine	57	10.1 Packaging	69
6 Application and Operation	59	10.2 Waste equipment	69
6.1 Switching between wet and dry processing	59	10.3 Coolant/Lubricant (C/L)	69
6.2 Inserting the blank holders into the blank changer (for Ceramill Matron only)	60	11 Technical Data, Spare Parts and Accessories	70
6.3 Changing the blank holders	60		
6.4 Program execution	61		
6.5 Checking the blank during a running program	61		
6.6 Prioritizing a newly inserted holder during a milling process (for Ceramill Matron only)	61		

1 Explanation of Symbols

Warning indications

 Warning indications in the text are marked with a triangle and boxed.

 In case of hazards through electricity, the exclamation mark in the warning triangle is substituted by a lightning bolt.

Signal words at the beginning of a warning indication specify the type and severity of the consequences, if the measures to avert the hazard are not adhered to.

- _ **NOTE** means that property damage can occur.
- _ **CAUTION** means that light to fairly serious personal injury can occur.
- _ **WARNING** means that serious personal injury can occur.
- _ **DANGER** means that life-threatening personal injury can occur.

Important information

Important information that do not lead to hazards for humans or property damage are marked with the icon aside and are also boxed.

Other symbols in the Manual

Symbol	Meaning
>	Item of an operation description
—	Item of a list
▪	Subitem of an operation description or a list
[3]	Numbers in square brackets refer to position numbers in graphics/figures

Other symbols on the machine

Symbol	Meaning
	USB port
	Network connection (Ethernet)
	Control cable connection, extraction
	Fuse
	Compressed-air connection
	Control cable connection, tank
	Supply connection, coolant/lubricant
	Connection for cleaning gun

EN

2 Safety Instructions

2.1 General Safety Instructions

When using the product, always observe the following safety instructions:

DANGER:

Danger of electric shock from leaking coolant/lubricant.

- ▷ When changing the coolant/lubricant, use the intended pumping device.
- ▷ If coolant/lubricant has leaked out: Disconnect the machine from the power supply and contact the technical service.

WARNING:

Health hazard during dry processing without extraction. Without suitable extraction, hazardous respirable milling dust can escape.

- ▷ Only carry out dry processing with suitable extraction.

CAUTION:

Damage of slipping from leaking coolant/lubricant!

- ▷ Regularly check the coolant/lubricant connection for tightness against leaks.
- ▷ When changing the coolant/lubricant, use the intended pumping device.
- ▷ If coolant/lubricant has leaked out: Set up warning signs and clean floor properly.

CAUTION:

Danger of tipping over when the drawer for the tools and tank is pulled out.

- ▷ Do not use the drawer as a climbing aid or seat.

CAUTION:

A defective product can lead to malfunctions!

On detection of damage or a functional defect of the product:

- ▷ Label the product as defective.
- ▷ Prevent further operation until the product has been repaired.

NOTE:

The system may only be used for the purpose that it is intended for. When disregarding, damage to property and personal injury cannot be excluded. Any and all warranty, guarantee and responsibility on behalf of the manufacturer shall expire.

NOTE:

The system may only be operated in enclosed rooms.

**NOTE:**

Possible damage due to insufficient load-bearing capacity of the floor.

The set-up location and the surrounding paths must be able to withstand the weight of the system.

- ▷ Observe weight of system.
- ▷ Inform yourself about the building conditions at the responsible authority (e.g. architect/structural engineer).

**NOTE:**

Possible machine damage through insufficient extraction during dry processing!

- ▷ Operate the machine only with the extraction system recommended by the manufacturer or another compatible extraction system.

**NOTE:**

Possible machine damage through insufficient cooling during wet processing!

- ▷ For wet processing, operate the machine only with the tank specifically intended for it.

**NOTE:**

Possible machine damage due to insufficient or incorrect compressed-air supply!

- ▷ Operate the machine with the required pressure, air volume and air quality (see Technical Data).

**NOTE:**

Switch the machine off when not in use or unsupervised for longer periods, e.g., overnight. This method of saving electrical energy also benefits the environment.

2.2 Unattended operation**NOTE:**

The machine may be operated unattended, provided that national and local laws and regulations permit this and these are complied with. Exception: Processing of titanium blanks (see section 2.3). Furthermore, the requirements of the respective insurance company must be met. The following conditions must be observed to ensure unattended operation:

- ▷ Do not use the machine with a heavily soiled milling chamber.
- ▷ Protect the machine against unauthorized access.
- ▷ Ensure that the room in which the machine is operated has an automatic fire alarm system.

2.3 Processing titanium blanks**NOTE:**

Despite water cooling, processing titanium involves an increased risk of fire. Metal shavings/chips can easily ignite and cause a metal fire.

For titanium indications:

- ▷ Do not operate the machine unattended.
- ▷ Have suitable extinguishing agents readily available.

3 Suitable Personnel

NOTE:

The machine may only be unpacked, installed or relocated by trained specialist personnel, authorized by the manufacturer.

NOTE:

The machine may only be put into operation by trained specialist personnel, authorized by the manufacturer.

NOTE:

The machine may only be operated by trained specialised personnel.

- _ Torque wrench, 2 Nm
- _ Blank holder 98 mm
- _ Hand scanner
 - Scanner
 - USB radio receiver
 - Charging cable
- _ Extractor adapter
- _ Extraction socket
- _ Plug for dry processing
- _ Cleaning brush, small
- _ Cleaning brush, large
- _ Terry towel

Optional scope of delivery:

- _ Torque wrench for collet
- _ Collet, 3 mm
- _ Auto-calibration kit
 - Auto-calibration holder
 - Calibration pin
 - Spiral cable

4 Product Specifications

4.1 Processing machine

This User Manual describes different models of the processing machine. Therefore, the representations can deviate from the actual appearance of the machine.

4.1.1 Delivery Scope

- _ Digital processing machine
- _ Mains cable
- _ USB cable (3 m)
- _ Network cable (5 m)
- _ Service unit incl. compressed-air hoses
- _ Service set for spindle
- _ Set for manual calibration
 - Test blank
 - Cutter Roto RFID grey 1.0
 - Cutter Roto RFID Calibration EM
 - Slide calipers

To check the zero point, a test-milling run was carried out at the factory before delivery. Milling chips may have remained from this. These are not a reason for a complaint, as this inspection process contributes significantly to the precision of the machine.

- ▷ After unpacking, check the machine for completeness and possible transport damages. Please claim any transport damages immediately with your supplier.
- ▷ After installation by the authorised service technician: Check scope of delivery and correct function. In case of missing parts or incorrect function, immediately file complaint on site with the service technician.

4.1.2 Intended Use

The processing machine is a CNC-controlled milling and grinding machine for the fabrication of dentures by means of:

- _ Dry processing of blanks made of
 - pre-sintered zirconium oxide
 - unsintered non-precious metals (NPMs)
 - Cobalt-chromium alloys
 - Wax
 - Plastics
 - Modelling materials
- _ Wet processing of blanks made of
 - Sintered glass ceramics
 - Plastics
 - Hybrid materials
 - Feldspathic ceramics
 - Titanium
 - Cobalt-chromium alloys

Using blanks and tools not approved by the manufacturer can damage the machine and make the product unusable. For such cases, the manufacturer shall assume no liability whatsoever.

Unauthorised modifications/alterations of the machine shall also void the warranty.

A compressed-air connection is required for operation of the processing machine (see page 54).

Dry processing is permitted only in conjunction with an appropriate extraction system (for technical requirements, see page 54).

Wet milling and grinding is permitted only in conjunction with Ceramill Coolantflow (accessory).

4.1.3 CE Declaration of Conformity

This product was designed and manufactured based on careful selection of the harmonised standards to be observed, as well as additional technical specifications. It thus corresponds with the state-of-the-art and ensures maximum safety.

In terms of design and performance, this product complies with the European Directives and the supplementary national requirements. Conformity has been confirmed with the CE marking.

Directives:

- _ 2006/42/EC (Machinery Directive)
- _ 2011/65/EU (RoHS Directive)
- _ 2004/108/EC (EMC Directive)

Harmonised standards:

- _ DIN EN IEC 63000:2019
- _ UL 61010-1:2012
- _ CAN/CSA-C22.2 Nr. 61010-1:2012
- _ DIN EN 61010-1:2020
- _ ISO 16090-1:2022

The product's declaration of conformity is available upon request from the manufacturer.

4.1.4 Components and interfaces

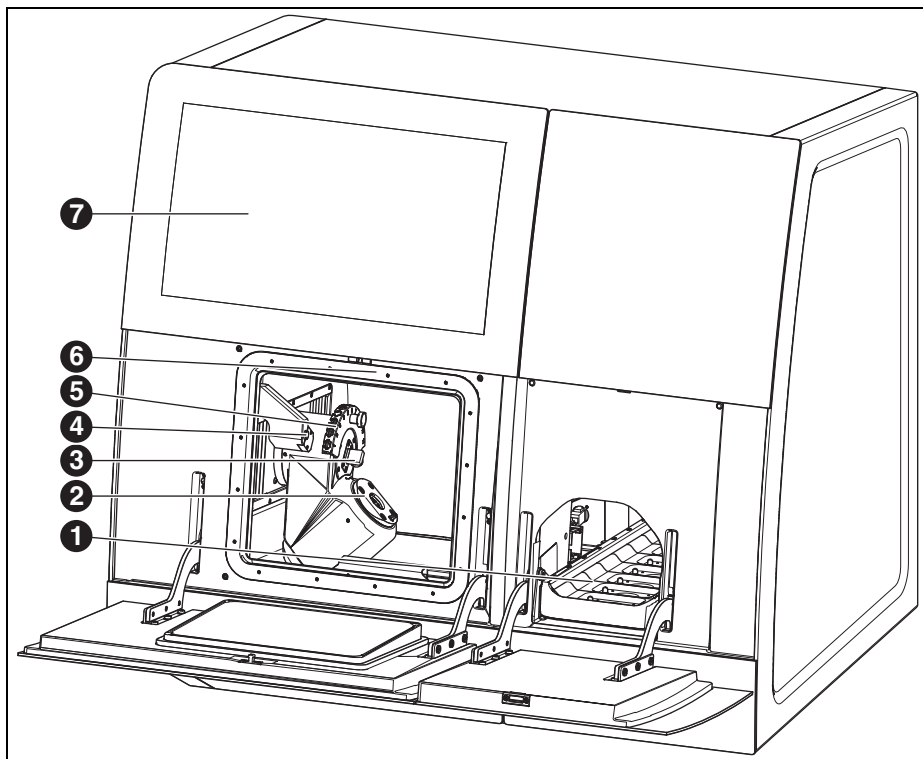
Fig. 1 Ceramill Matron front

- 1 Blanks magazine
- 2 Zero-point clamping system
- 3 Measuring device, tool length
- 4 Presence sensor/RFID antenna
- 5 Magnetic tool holder (SnapMag)
- 6 Covered: Milling spindle
- 7 Touch display

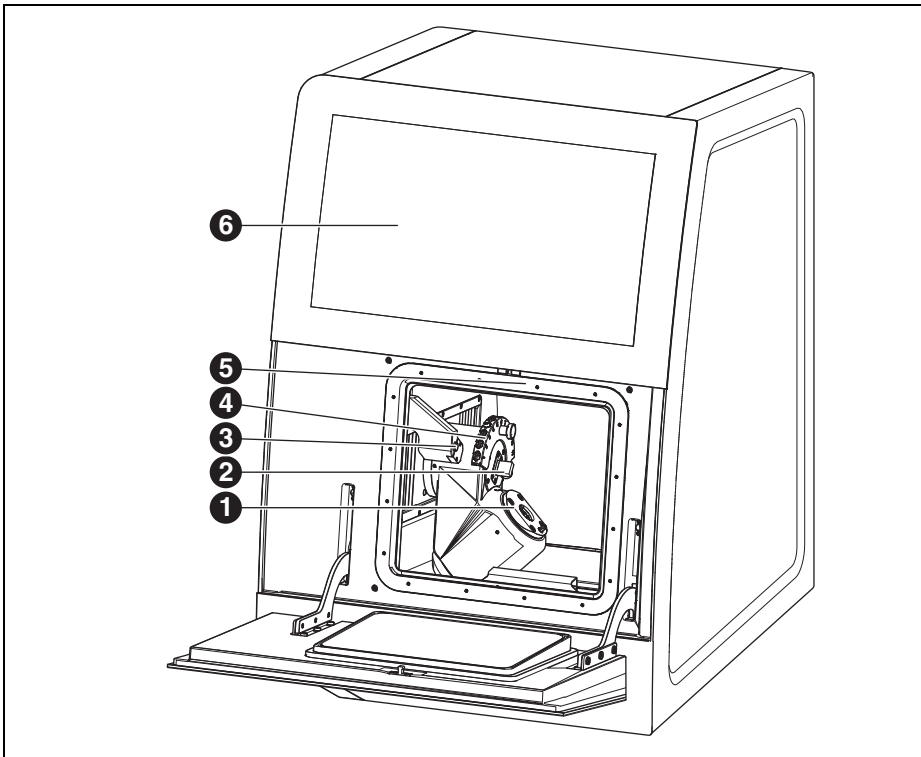
Fig. 2 Ceramill Matron Core front

- 1 Zero-point clamping system
- 2 Measuring device, tool length
- 3 Presence sensor/RFID antenna
- 4 Magnetic tool holder (SnapMag)
- 5 Covered: Milling spindle
- 6 Touch display

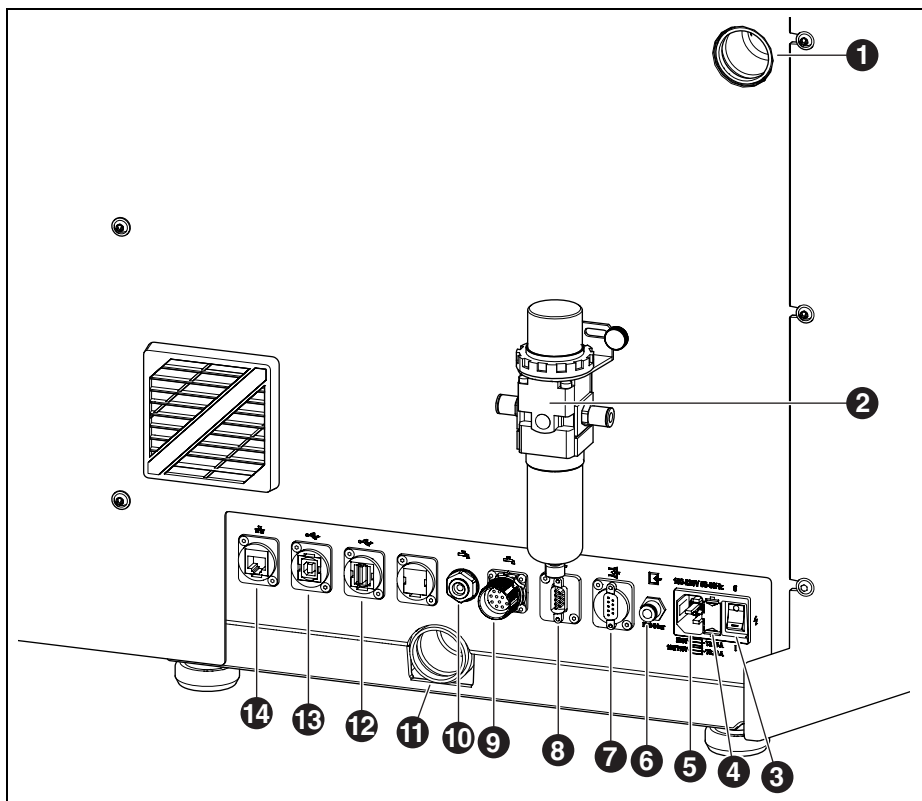
Fig. 3 Ceramill Matron/Ceramill Matron Core rear

- | | |
|--|---|
| 1 Extraction-hose connection | 13 USB-B port  |
| 2 Service unit | 14 Ethernet interface  |
| 3 Main switch | |
| 4 Drawer with fuse holder | |
| 5 Power supply connection | |
| 6 Compressed-air connection (of service unit) | |
| 7 Control interface for extraction  | |
| 8 Control interface for extraction  | |
| 9 Control interface for pump unit  | |
| 10 Connection, C/L supply line  | |
| 11 Connection, C/L drain | |
| 12 USB-A port  | |

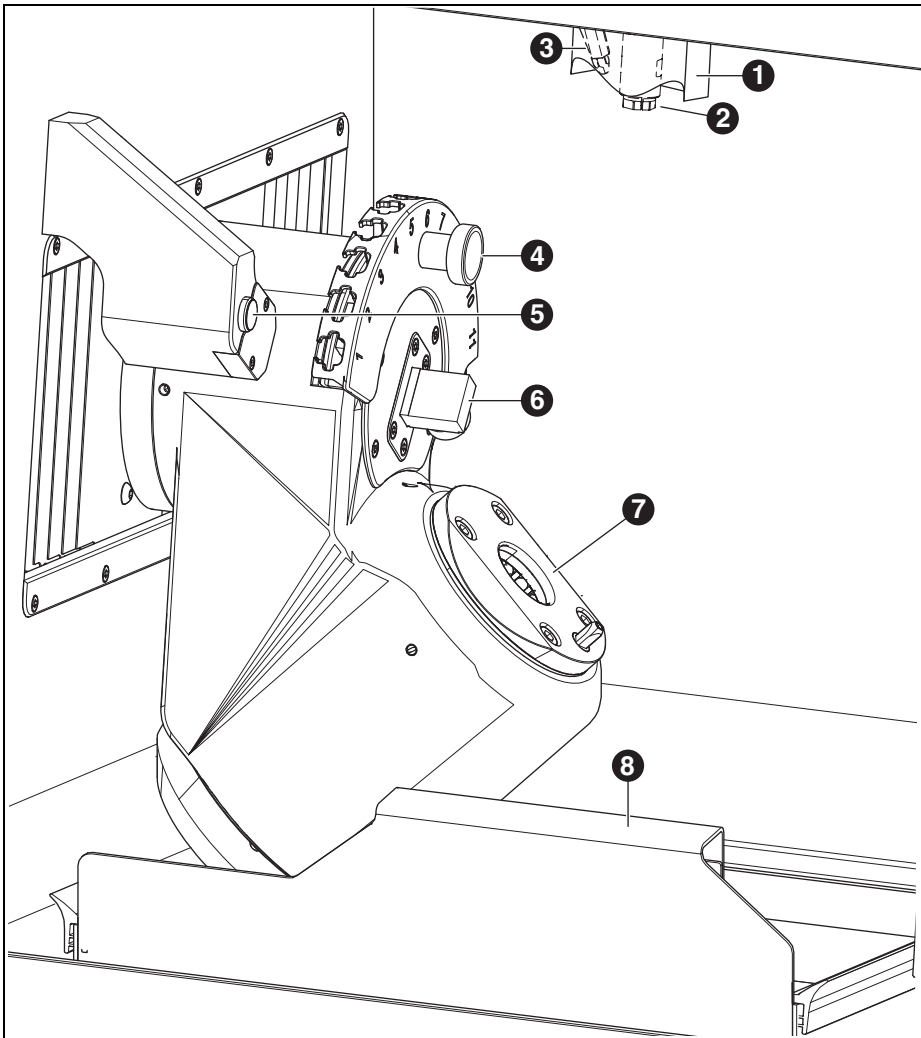
Fig. 4 Milling chamber

- | | |
|--|---------------------------------|
| 1 Extraction socket | 6 Measuring device, tool length |
| 2 Spindle | 7 Zero-point clamping system |
| 3 Compressed-air nozzle (for dry processing) | 8 Chip tray with sieve |
| 4 Magnetic tool holder (SnapMag) | |
| 5 Presence sensor/RFID antenna | |

4.2 Accessories

4.2.1 Milling and grinding tools

Various cutting and grinding tools are available for the processing machine (see catalogue).

4.2.2 Blank holder

Various blank holders are available for the processing machine (see catalogue).

4.2.3 Extraction

An extraction recommended by the manufacturer is required for dry processing with the processing machine (for technical requirements, see page 70).

4.2.4 Tank

Intended Use

The tank enables wet processing in the processing machine. It ensures that the required coolant/lubricant is drained and filtered in accordance with the application.

Operation with machines other than the processing machine is not permitted.

Delivery Scope

- _ Tank with
 - Container for coolant/lubricant (C/L)
 - Filter for (C/L)
- _ C/L supply line
- _ Piping kit
- _ 2 litres of coolant/lubricant
- _ Analytics set
 - Refractometer
 - Test strips for pH value
 - Test strips for nitrite content
- _ Instructions

Components and interfaces

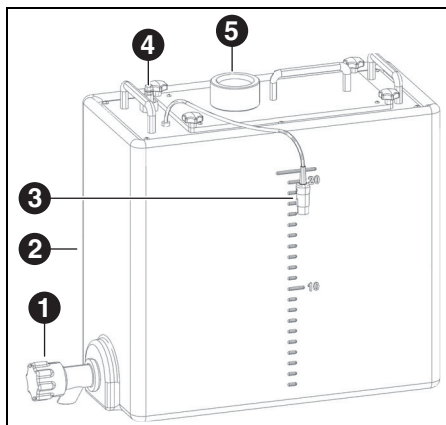
Fig. 5 Tank connections

- 1 C/L drain cock
- 2 C/L container
- 3 Control interface for the tank
- 4 Connection to pump unit
- 5 Feed from the processing machine

4.2.5 Pump unit (CoolantFlow)

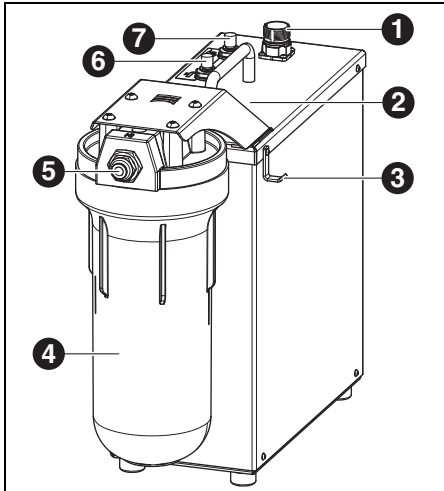
Fig. 6

- 1 Control interface with the processing machine 
- 2 Control interface for tank (covered)
- 3 Hook for cleaning gun
- 4 Filter
- 5 C/L inlet
- 6 C/L outlet to processing machine 
- 7 C/L outlet to cleaning gun 

Intended Use

In combination with the tank, the pump unit enables wet processing in the processing machine. It ensures the application-efficient supply, drainage and filtering of the required coolant/lubricant.

Operation with machines other than the processing machine is not permitted.

Delivery Scope

- _ Pump unit with
 - Filter set
 - Cleaning gun with hose
 - Chip tray
 - Connection cable
- _ Connection hose from tank to pump unit
- _ Connection hose from pump unit to machine
- _ Drain hose from machine to tank

EN

4.2.6 Cart

The cart serves as a base frame for the processing machine. The feet of the processing machine are positioned into the retainers of the cart. The tank, pump unit and extraction system can be installed in the base assembly of the cart.

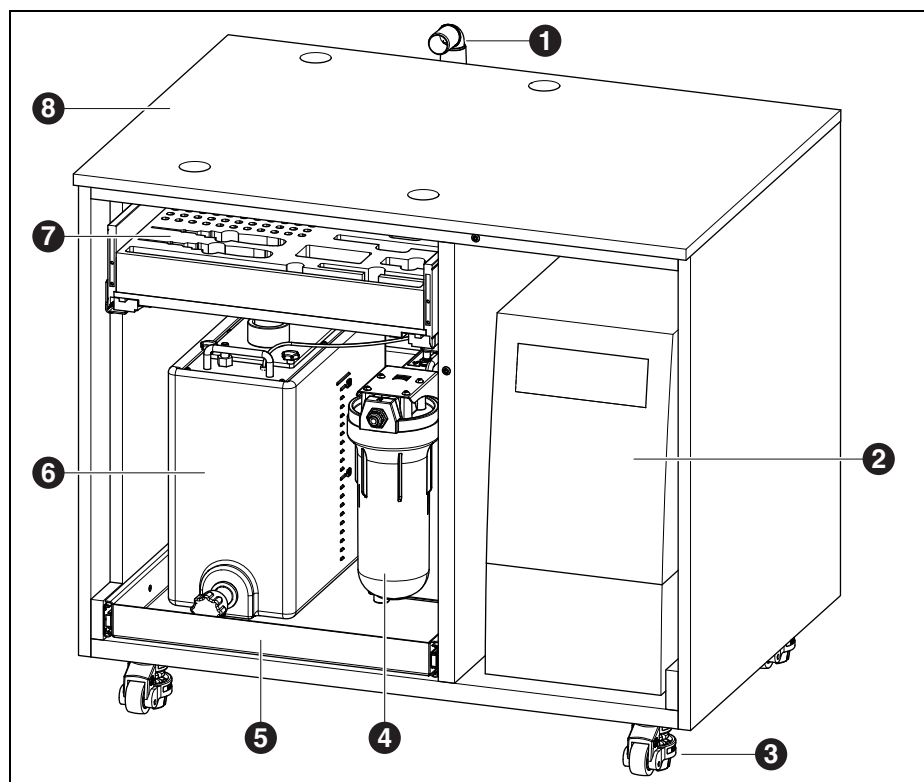
Fig. 7 Cart (doors hidden)

- | | |
|---------------------------------|------------------------------------|
| 1 C/L drain | 7 Drawer for tools and accessories |
| 2 Extraction | 8 Cart |
| 3 Castors, lockable (4 x) | |
| 4 Pump unit | |
| 5 Drawer for tank and pump unit | |
| 6 Tank | |

5 Installation

5.1 Processing machine

Countries with 100 - 115 V mains voltage

The machine is preset for a mains voltage of 230 V. When the mains voltage at the set-up location is 100 - 115 V:

- ▷ Pull out the drawer ([4] in Fig. 3 on page 46) next to the mains plug.
- ▷ Replace both T3.15 A fuses with the provided T6.3 A fuses.

5.1.1 Setting up the processing machine

Please refer to the separate processing machine set-up instructions.

- _ The machine is intended exclusively for use within dry, closed rooms.
- _ A clearance of at least 200 mm is to be observed sideways, to the rear and upwards. Minimum space requirements for the processing machine (incl. connections):

	Ceramill Matron	Ceramill Matron Core
Width	1340 mm	974 mm
Depth	850 mm	850 mm
Height	1050 mm	1050 mm

Tab. 1

For maintenance, 1 meter of free space is required in front of, behind and next to the machine. This free space only needs to be present on one side if the other sides can be made accessible by turning the machine.

- _ The empty weight of the Ceramill Matron is 150 kg, and that of the Ceramill Matron Core is 110 kg. In addition, the weight of the coolant/lubricant in the container and the blanks including the holder in the changer must be taken into account, depending on the application. The set-up surface must withstand an appropriate load.
- _ During operation, the room temperature must be between 18 °C and 30 °C; avoid large variations in temperature.
- ▷ Set up the machine on a sturdy work table or work bench (off the ground and clear of walls), or on the cart available from the manufacturer.
- ▷ On the Ceramill Matron: Turn the support screws of the changer unit downwards to secure the processing machine against tipping over.
- ▷ Mount the service unit (see page 53).
- ▷ Insert the extraction hose of the extraction system to the connection on the rear side of the processing machine.

Setting up the processing machine on the cart

A clearance of at least 200 mm is to be observed sideways, to the rear and upwards. Space requirements for the processing machine on the cart incl. connections:

- _ Height: 1885 mm
- _ Width: 1340 mm
- _ Depth: 1150 mm

The set-up surface must be level.

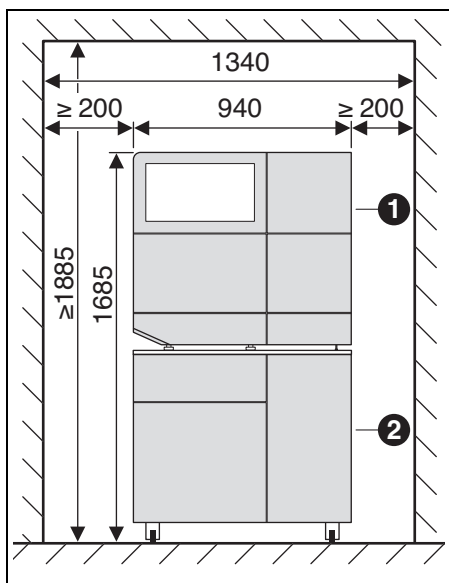
Fig. 8 Ceramill Matron on cart

- 1 Processing machine
- 2 Cart

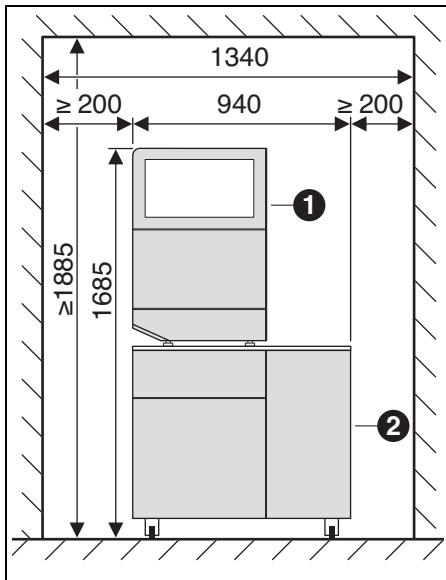
Fig. 9 Ceramill Matron Core on cart

- 1 Processing machine
- 2 Cart

NOTE:

Possible damage of the electric lines and the C/L supply line:

- ▷ When moving the cart, pay attention that no lines are driven over.

- ▷ On the Ceramill Matron: Turn the support screws of the changer unit downwards to secure the processing machine against tipping over.
- ▷ Position the processing machine onto the cart, ensuring that the legs engage into the retainers.
- ▷ Move the cart with the processing machine to the set-up location.
- ▷ Rotate the 5 feet down until they have firm contact.

5.2 Setting up the tank

Place the tank under the set-up table of the processing machine or into the processing machine cart, which is available from the manufacturer.

In order for the C/L drain pipe to be mounted, the table must be set up with a certain clearance to the wall (see set-up instructions).

- ▷ Fill coolant/lubricant into the C/L container (see section 7.5 on page 64).

NOTE:

Possible damage of the electric lines and the C/L supply line:

- ▷ When moving the cart, pay attention that no lines are driven over.

- ▷ Connect the C/L supply line to the processing machine and pump unit.
- ▷ Attach the C/L drain pipe and connect it to the processing machine and tank (see installation instructions).
- ▷ Connect the control cable to the control interface for the pump unit on the rear side of the processing machine ([9] in Fig. 3) and to the pump unit ([1] in Fig. 6).
- ▷ Set up the tank or install it into the cart (see section 7.6) as required.
- ▷ Connect interface ([3] in Fig. 5) with pump unit ([2] in Fig. 6).
- ▷ Connect the connection to pump unit ([4] in Fig. 5) to the C/L inlet ([5] in Fig. 6).

5.3 Service unit

The spindle is equipped with a sealing-air feature. This sealing air prevents shavings/chips and dust from entering the spindle. The service unit regulates the inlet pressure to the factory-set pressure. The integrated filter filters particles to a size of 5 µm.

NOTE:

The service unit does not protect against spindle damage caused by moist or contaminated compressed air.

- ▷ Ensure that the compressed air meets the required quality.

Installing the service unit

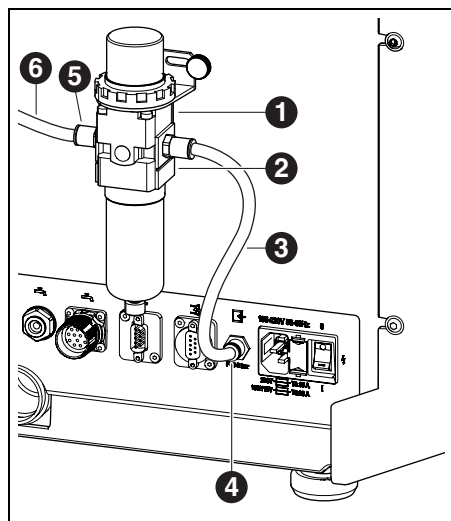
Fig. 10 Service unit

- 1 Service unit
 - 2 Pressure reducer
 - 3 Compressed-air hose, 250 mm
 - 4 Compressed-air connection
 - 5 Compressed air inlet of service unit
 - 6 Compressed-air hose, 2000 mm
- ▷ Mount the service unit on the rear side of the machine with the provided screws.
 - ▷ Connect the 250 mm compressed-air hose [3] to the right-hand connection of the service unit and to the compressed-air connection [4] of the processing machine.
 - ▷ Connect the 2000 mm compressed-air hose [6] to the compressed-air inlet [5] of the service unit and the other end to the compressed-air supply of the set-up location.

Requirements for the compressed-air supply

For operation, the compressed air must fulfil the following minimum requirements:

- _ Flow rate: 50 l/min (1.77 CFM)
- _ Minimum pressure: 6 bar/87 psi

Requirements on air purity:

- _ Solid contaminants: class 3;
 - Maximum number of particles per cubic meter of compressed air:
 - 90000 particles in the range of 0.5 – 1 μm
 - 1000 particles in the range of 1 – 5 μm
- _ Water content: class 4;
 - Max. pressure dew point +3 $^{\circ}\text{C}$
- _ Total oil content: class 2;
 - Max. oil content 0.1 mg/m^3

In continuous operation, the pressure is monitored by the processing machine. When the required values are fallen below, the milling program stops. Once the correct air values are available again, the milling program continues automatically.

Recommended extraction system requirements

If an extraction system not recommended by the manufacturer is used, this extraction system must fulfil the following requirements:

- _ Extraction capacity: 240 m^3/h
- _ Suitability for zirconium-oxide dusts
- _ HEPA microfilter (captures 99.97 % at a particle size of 0.3 μm)


The connection of the extraction hose to the machine is made with an adapter (in delivery scope). The adapter fits extraction hoses with a 38 mm interior diameter.

5.4 Installation of the CAM and machine software

An AG.LIVE account is required to receive and install Ceramill software products.


- For new customers: Register with a lab account in AG.LIVE.
- In the AG.LIVE dashboard, click on *Check for updates*.

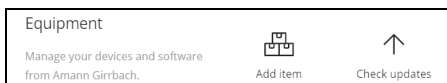The software products are indicated.



- Download and install the *Ceramill Software Manager*.
 - Open the *Ceramill Software Manager* and log in with the AG.LIVE access data.
 - In the Ceramill Software Manager open the *Available software* area.
- All available Ceramill software products are displayed.

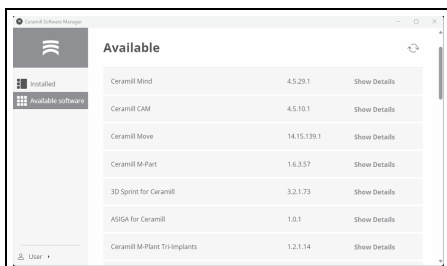
- Install the following Ceramill software products to operate the Ceramill Matron:

- Ceramill CAM
- Ceramill Move

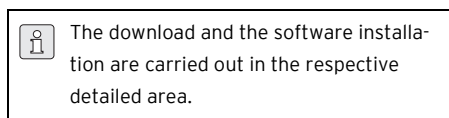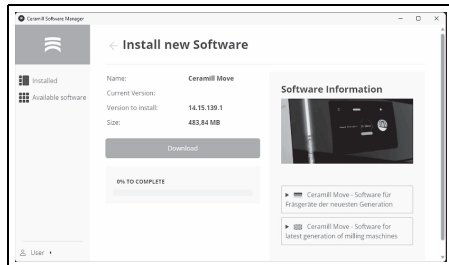After successful installation, the software products are in the *Installed* area.

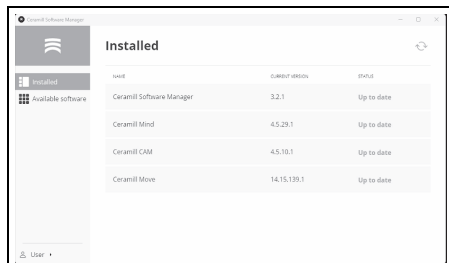5.5 Connecting the interfaces

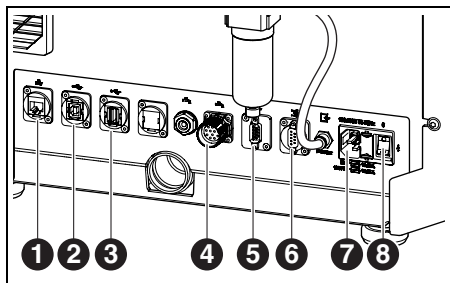
Fig. 11 Interfaces on the rear side of the processing machine

- 1 Ethernet interface
- 2 USB-B port
- 3 USB-A port
- 4 Control interface for pump unit
- 5 Control interface for extraction
- 6 Control interface for extraction
- 7 Power supply connection
- 8 Main switch

- ▷ Make sure that the main switch [8] of the processing machine is in the 0 position.
- ▷ Connect the mains cable to the power supply connection [7] of the processing machine and to a power outlet.
- ▷ Connect the control cable to the control interface for extraction [5] respectively [6] on the rear side of the processing machine and the rear side of the extraction system.

For installation with pump unit:

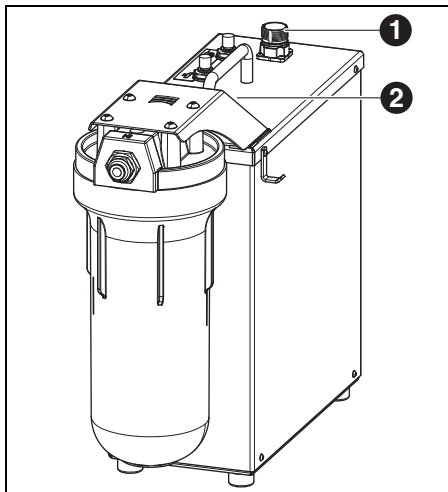
Fig. 12

- 1 Control interface with the processing machine
- 2 Control interface for tank (covered)

- ▷ Connect the control cable to the control interface for the pump unit on the rear side of the processing machine ([4] in Fig. 11) and to the pump unit ([1] in Fig. 12).
- ▷ Connect interface at tank ([3] in Fig. 5) with pump unit ([2] in Fig. 12).

5.6 Connection to the PC

- ▷ Connect a USB cable to the USB-B port ([2] in fig. 11) of the processing machine and to a free USB port of the PC.
- The maximum permitted length of the USB cable is 5 m.

5.7 Establishing the network connection

5.7.1 Ethernet

- Connect the Ethernet interface ([1] in Fig. 11) of the processing machine with the PC or with the network using a network cable.
The maximum permitted length of the network cable to the PC or to the network router is 100 m.
- Switch the processing machine on.
- Start the processing machine software.
The machine automatically recognizes the connection through which it is connected to the network.

The PC and the processing machine must both be in the same sub-network.

5.7.2 WLAN

- Plug the WLAN USB antenna into the intended slot on the rear side of the machine.
- Switch the processing machine on.
- Select the WLAN network on the display under *Settings*.
- Start the processing machine software.
The machine automatically recognizes the connection through which it is connected to the network.

5.8 Calibration of the processing machine

A calibration of the processing machine is mandatory after installation for commissioning. The processing machine can adjust itself via auto-calibration (optional). This process can be started on the display under Maintenance. In this, all necessary steps are explained on the display. Next up is a description of the calibration via the milling process. This must be carried out monthly or on request of the support.

The tools can be inserted in any of the SnapMag slots. The integrated RFID chip enables the machine to automatically recognize the tool. When a tool has reached its wear limit, the processing machine automatically accesses the respective sister tool (as far as available).

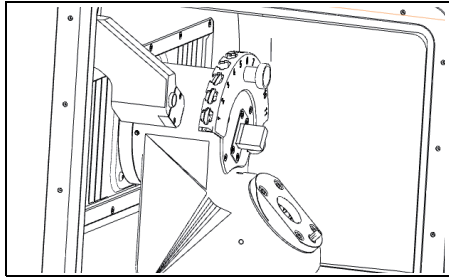
Fig. 13 Tool holder

- For calibration, insert the required cutters from the delivery scope into the tool holder.

- ▷ Insert the required test blank into the size 98 blank holder.

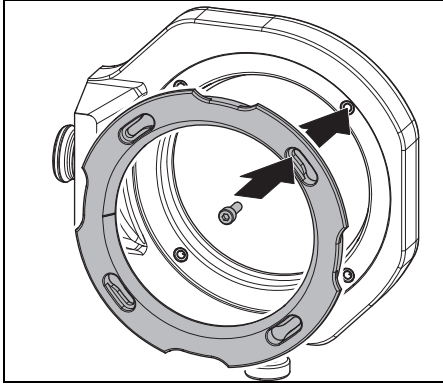
Fig. 14 Inserting a blank

- ▷ Fasten the blank using four screws.

NOTE:

Do not damage the blank!

- ▷ Do not overtighten the screws! Avoid mechanical stress!

- ▷ Tighten the screws with the enclosed torque wrench.

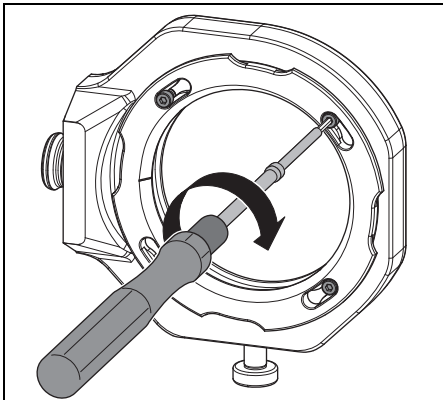
Fig. 15 Tightening a blank

The blank holder must either be in the zero-point clamping system or in the blank changer.

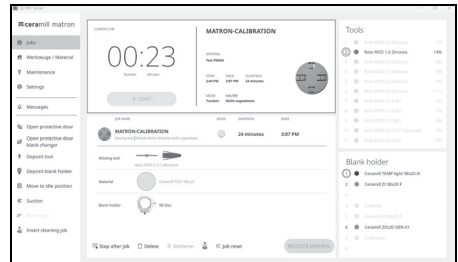
Fig. 16 Screenshot Linking the holder

- ▷ Link the holder to the blank at the display.
- ▷ Start the calibration process on the display under *Maintenance*.

The *Calibration* workflow guides you through all necessary steps for a correctly set processing machine.

6 Application and Operation

WARNING:

Danger from shavings/chips or tool fragments being thrown from the machine!

- ▷ Make sure that the door of the processing machine is always closed during processing!

NOTE:

Insufficient milling results and increased wear when dry processing!

- ▷ Do not operate the processing machine without extraction!

NOTE:

Insufficient milling/grinding results and increased wear when wet processing!

- ▷ Do not operate the processing machine without supply unit!

6.1 Switching between wet and dry processing

The processing machine can be used both for wet and dry processing.

In delivery condition, it is set up for dry processing.

6.1.1 Switching from dry to wet processing

- ▷ Remove chip tray with sieve and clean.
- ▷ Remove extraction socket ([1] in Fig. 4) from the extraction opening of the spindle.
- ▷ Connect the extraction hose to the extraction opening.
- ▷ Activate the extraction system on the display.
- ▷ Vacuum the milling chamber.
- ▷ Clean the tool holder.

When the milling chamber is clean:

- ▷ Deactivate the extraction system on the display.
- ▷ Remove the extraction hose.
- ▷ Attach extraction nozzle to spindle.
- ▷ Insert the cleaned chip tray with the sieve facing to the rear.

For wet processing, it is mandatory to use the chip tray. The sieve serves as the first coolant/lubricant filter stage for coarse particles.

- ▷ Check if the nozzles for the coolant/lubricant are free of contamination.
- ▷ Check if the door gasket is clean and free of swarf/chips.
- ▷ Check if the C/L drain pipe is free of obstructions and if the piping on the outside is securely connected.
- ▷ Start the milling process.

6.1.2 Switching from wet to dry processing

- ▷ Remove the chip tray.
- ▷ Clean and dry the milling chamber and chip tray (see page 59).
- ▷ Clean the tool holder.

The chip collection tray [3] is not necessary for dry processing, but it can remain in the milling chamber.

- ▷ Check the extraction.
- ▷ Start the milling process.

6.2 Inserting the blank holders into the blank changer (for Ceramill Matron only)

In order for the processing machine to be able to change the blank holder automatically, it must be placed correctly in the blank changer.

- ▷ Open the door of the blank changer via the menu on the display.
- ▷ Pull out the blank-holder tray with the blank-holder retainer.
- ▷ Insert the blank holder with the RFID chip to the right (from the front).
- ▷ Push the blank -holder tray back in until beyond the latching point (can be felt).
- ▷ Close the door.

The magazine is automatically lowered. Once all axes are at rest, an inventory run begins.

The blank holder must not become wedged in the retainer.
When the tray has not been pushed over the latching point, the processing machine outputs an error indication.

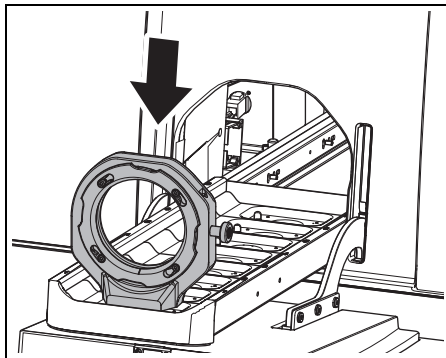
Fig. 17 Inserting a blank holder into the blank changer

6.3 Changing the blank holders

In delivery condition, the processing machine is equipped with a size 98 blank holder. For processing varying materials, it may be required to use a different blank holder.

On the Ceramill Matron:

- ▷ Affix blank in blank holder and insert it into the blank changer.

The blank holders are changed automatically via the blank changer after the blank has been assigned to the holder and the job.

On the Ceramill Matron Core:

- ▷ Affix blank in blank holder.
- ▷ Press *Change blank holder* on the machine display.
- ▷ Follow the instructions on the display.

6.4 Program execution

Before executing a program, it must be transferred to the machine (see CAM software). After transferring, a new job is indicated on the display of the processing machine.

- ▷ On the display of the processing machine, start the new job with Start job.

Several jobs can be transferred at the same time. These are then arranged in the order in which they were selected on the PC in a job list.

The program can be stopped:

- ▷ On the display of the processing machine, press Pause job.

To continue the program:

- ▷ With the milling chamber door closed, press Continue job on the display.

The program is continued at the position where it was stopped.

At the end of the program, the machine performs an automatic cleaning and then goes into the stand-by position or starts the next job.

6.5 Checking the blank during a running program

- ▷ On the display of the processing machine, press Pause job.

The program is interrupted. The spindle moves to the rear.

- ▷ Open the milling chamber door.

- ▷ Visually check the blank.

- ▷ Close the milling chamber door.

- ▷ On the display of the processing machine, press Continue job.

The tool is picked up again. The spindle travels to the position, where the program was interrupted.

6.6 Prioritizing a newly inserted holder during a milling process (for Ceramill Matron only)

- ▷ Insert the blank holder into the blank changer and prioritize the job via the display.

- ▷ On the display of the processing machine, press Pause job.

The program is interrupted. The blank changer carries out an inventory.

- ▷ After completion of the inventory, press Continue job on the display of the processing machine.

The tool is picked up again. The spindle travels to the position, where the program was interrupted and continues the interrupted milling process.

7 Handling the Coolant/Lubricant (C/L)

CAUTION:

Health hazards are possible through improper handling of C/L!

- ▷ Use only C/L recommended by the manufacturer!
- ▷ Observe the safety instructions and handling regulations for C/L!

7.1 C/L analysis

For legal and technical reasons, the coolant/lubricant must be analysed weekly.

- ▷ Observe country-specific laws and guidelines on coolant/lubricant!

Regular analysis of the C/L is required for proper function of the processing machine. The results of the analyses must be documented and archived.

- ▷ Check the following coolant/lubricant properties weekly:
 - Concentration
 - pH-value
 - Nitrite/nitrate content

The analysis is to be carried out with test strips and refractometer. The test strips are simple to handle and enable the pH-value and the nitrite/nitrate content to be determined. The concentration of the C/L is measured with the refractometer.

- ▷ On the display, call up the Check coolant workflow.
- ▷ Place a clean, flat receptacle in the centre of the milling chamber.
- ▷ Use the cleaning gun to fill the receptacle to a height of approx. 2 cm.

- ▷ Remove the receptacle.
- ▷ Measure the pH and nitrite/nitrate values:
 - Immerse the reaction zone of the pH and nitrite test strip into the C/L.
 - Shake off excess fluid.
 - After one minute, compare the reaction zone with the colour scale and read off the measured value.
- ▷ Measuring the concentration:
 - Use the liquid dropper to remove the C/L and apply it to the measuring surface of the refractometer (glass zone).
 - Close the refractometer lid.

The measuring surface must be completely wetted with liquid.
 - Hold the refractometer against a light source and look through the ocular (eyepiece).

The reference values and the information on measures in case of exceeding or falling below the reference values can be found on an information sheet in the delivery scope of the tank.

7.2 Removing the tank from inside the cart

NOTE:

Possible damage of the electric lines and the C/L supply line:

▷ When moving the cart, pay attention that no lines are driven over.

- ▷ Open the cart door and pull out the drawer with the tank and pump unit.
- ▷ Pull the C/L drain pipe [1] and the C/L supply line [3] upward off of the tank.
- ▷ Pull off control cable [2] from pump unit.

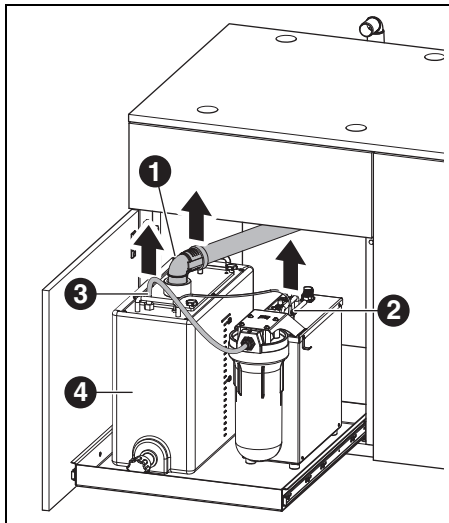
Fig. 18 Removing the tank interfaces

- 1 C/L drain pipe
- 2 Control cable
- 3 C/L supply line

▷ Carefully lift tank [4] out of the cart.

7.3 Emptying the C/L container

NOTE:

Possible damage to the drain cock!

▷ Move/carry the tank only by the handles!

▷ Set tank [1] onto a table with the drain cock facing toward the front.

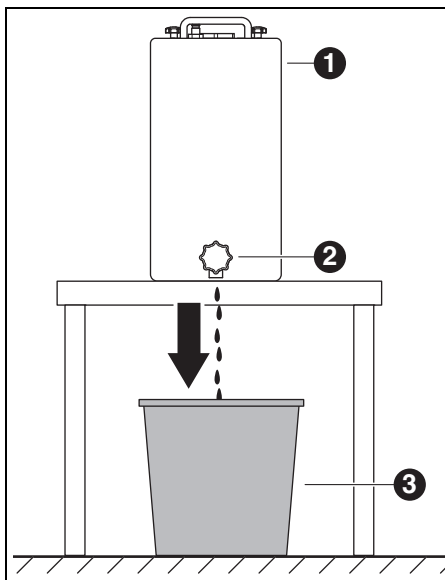
Fig. 19

- 1 Tank
- 2 Drain cock
- 3 Disposal receptacle

- ▷ Position a disposal receptacle [3] below the drain cock [2].
- ▷ If required, mount a hose to the drain cock and conduct it into the disposal receptacle.
- ▷ Open drain cock.

The C/L flows out of the tank into the disposal receptacle.

When the tank is empty:

- ▷ Shut the drain cock again.

The coolant/lubricant must be disposed of in accordance with the applicable regulations and laws. Proper disposal of the C/L must be documented.

7.4 Cleaning the tank

NOTE:

The filter and the tank may only be cleaned with water.

- ▷ Do not use cleaning agents.

- ▷ Unscrew and remove the four star-knob nuts on the top side of the tank.
- ▷ Remove the lid with the drain pipe.
- ▷ Remove the filter bag and empty it.
- ▷ Clean the filter bag.
- ▷ Rinse out and clean the tank.

7.5 Filling in coolant/lubricant

The water for the coolant/lubricant must have drinking-water quality. The optimal hardness range is between 10 - 20 °dH (corresponds to 1.8 - 3.6 mmol total hardness per litre). If the hardness is far beyond this range, we recommend the addition of distilled water. For grinding glass-ceramic materials, it is generally recommended to use distilled water.

- ▷ Fill 19 litres of water and 1 litre of coolant concentrate into the tank and mix well.
- ▷ Insert filter bag into tank.
- ▷ Reattach the lid with the drain pipe and secure it with the four star-knob nuts.

7.6 Installing the tank into the cart

- ▷ Lift the tank into the cart.
- ▷ Connect the control cable to the pump unit.
- ▷ Connect the C/L drain pipe and the C/L supply line to the top to the tank.

NOTE:

Danger of damaging cables and hoses!

- ▷ When moving the cart, do not pinch or drive over any cables or hoses.

7.7 Cleaning the C/L system

- ▷ Remove chip tray and clean.
- ▷ Clean the milling chamber.
- ▷ Reinsert chip tray.
- ▷ Empty C/L container (see page 63).
- ▷ Clean C/L container (see page 64).
- ▷ Fill 5 litres of water and 250 ml of coolant concentrate into the C/L container and mix well.
- ▷ On the display, select Maintenance > Clean the C/L circuit.
- ▷ Rinse processing machine for 15 minutes.
- ▷ Empty the C/L container.
- ▷ Insert the filter bag.
- ▷ Fill 20 litres of water and 1 litre of coolant concentrate into the C/L container and mix well.
- ▷ Reconnect the tank completely.

7.8 Changing the filter cartridge

NOTE:

Possible blockages and damage to the C/L pump!

- ▷ Changing the filter cartridge regularly.

For regular wet processing, the filter cartridge of the pump unit must be replaced weekly.

- ▷ On the display, via Maintenance, call up the workflow for changing the filter cartridge and follow the instructions.

8 Cleaning and Maintenance

On the display, via Maintenance, interactive cleaning workflows can be called up.

The following applies for all cleaning workflows:

- _ Instructions are shown on the display when due.
- _ They are documented with pictures and can be called up and carried out step by step on the display of the processing machine.
- _ They are documented.

8.1 Processing machine

8.1.1 Cleaning

The processing machine must be cleaned after each working day.

- ▷ Empty the chip tray frequently, especially when wet milling plastic.

When the chip tray is too full, it is possible that the coolant/lubricant will not drain off.

Cleaning after dry processing

If there is residual zirconium oxide dust in the milling chamber from dry processing, this can clog the filter bag of the tank during subsequent wet processing.

NOTE:

Zirconia dust is abrasive!

- ▷ Clean the window in the door carefully, using only a soft cloth; do not scrub or scour.

- ▷ Open the door of the processing machine.
- ▷ Replace the extraction nozzle with the cleaning hose with brush head.
- ▷ On the display of the processing machine, activate Exhaustion.

- ▷ If necessary, remove the chip tray and dispose of the milling waste outside the machine.
This makes the extraction fill up less quickly, allowing for prolonged use.
- ▷ Thoroughly vacuum out the milling chamber with the extraction hose.
- ▷ Carefully clean the RFID read/write unit with a clean cloth.
- ▷ Clean tool shafts carefully.
- ▷ On the display of the processing machine, deactivate Exhaustion.

Cleaning after wet processing

NOTE:

Possible machine damage through dried-in coolant/lubricant and bound in milling dust.
After wet processing:

- ▷ Clean and dry the milling chamber prior to downtimes lasting several hours.

NOTE:

The extraction is not designed for wet vacuuming. The extraction may not be used for cleaning the interior after wet processing!

- ▷ Remove shavings/chips manually.

NOTE:

The milling chamber may only be cleaned with water.

- ▷ Do not use cleaning agents.

- ▷ Remove the chip tray and dispose of the milling waste in it.
- ▷ Clean the milling chamber and wipe it dry using a soft cloth.

Damage caused by not removed grinding residue are not covered by the warranty.

8.1.2 Regular inspections

- _ Condition of the C/L.
- _ Condition of hoses and drain pipes
- ▷ Keep the door gasket and the stopping surface on the door clean. Ensure not to damage the sealing lips.

8.1.3 Weekly maintenance

NOTE:

Danger of damaging the machine!

- ▷ Do not use compressed air, ultrasonic or a steam jet for the cleaning!

- ▷ On the display of the processing machine, select Maintenance > Clean machine.
- ▷ Follow the instructions on the display.

Spindle

For maintenance of the spindle ([2] in Fig. 4 on page 47), the enclosed spindle service set is required.

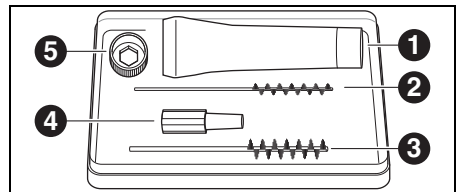
Fig. 20 Service set for spindle

- 1 Collet grease
- 2 Collet, 4 mm
- 3 Collet, 7 mm
- 4 Cone wiper
- 5 Collet key, size 10

In order to clean the spindle:

- ▷ On the display of the processing machine, select Maintenance > Clean chuck.

- ▷ Remove any dust from the spindle using a dry, clean cloth. Do not use cleaning agents!

Instead of the collet key, the collet torque wrench (optional accessory) can also be used for maintenance.

To remove the collet:

- ▷ Apply collet key [5] to collet.

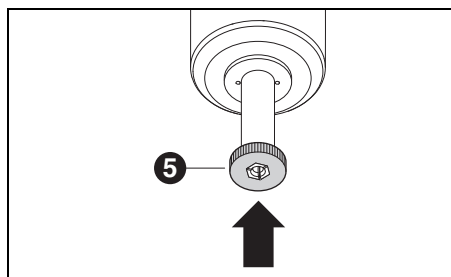
Fig. 21

- ▷ Using the collet key, unscrew the collet out of the spindle.

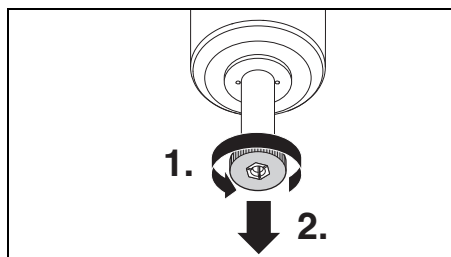
Fig. 22

- ▷ Clean collet holder with cone wiper [3].

The collet holder must be free of shavings/chips and debris.

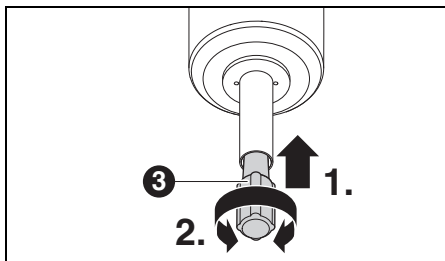
Fig. 23

- ▷ Clean the inside of the collet with the fine collet brush [2].

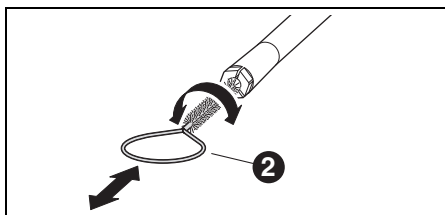
Fig. 24

- ▷ Apply a light coat of the collet grease ([1] in Fig. 20) from the service set to the outer cone of the collet.

This improves the sliding properties and increases the clamping force of the collet.

NOTE:

Only use the provided collet grease for lubrication of the collet cone. Do not use other types of grease.

- ▷ Screw the cleaned and lubricated collet in reverse order to the stop back into the spindle.

Service unit

- ▷ Visually check the service unit for:
 - Visible particles, debris or particulate matter within the condensate container (signs for contaminated compressed air)
 - Amber discolouration of the filter element (sign of oil in the compressed air)
 - Water in the filter (sign of water in the compressed air)

Once any of the listed impairments is given:

- ▷ Correct the problem with the compressed air.
- ▷ Replace complete service unit.

8.1.4 Monthly maintenance

Checking the zero point

We recommend carrying out a monthly zero point check using the optional "auto-calibration" accessory.

- ▷ On the display of the processing machine, select Maintenance > Calibrate machine and follow the instructions on the display.

8.1.5 External maintenance

External maintenance is necessary after max. 1000 operating hours. The processing machine will inform you appropriately.

8.2 Tank

8.2.1 Inspections and maintenance measures

Regular inspections:

- _ Condition of hoses and drain pipes
- _ C/L filling level in the tank.

When the filling level in the C/L tank has fallen below the minimum level, the milling program is interrupted and an error message is indicated.

- ▷ Fill up the tank (see page 64).
- ▷ Acknowledge error in the software.

8.2.2 Weekly maintenance

- ▷ Analysing the C/L (see section 7.1, page 62).

8.2.3 Half-yearly maintenance

- ▷ Clean the C/L system (see section 7.7, page 65).

The half-yearly cleaning is indicated by the program. Upon activation of the system cleaning, the indication is blanked out again.

9 Malfunctions, Repairs and Warranty

9.1 Malfunctions

In case of malfunctions:

- ▷ Restart the software.
- ▷ Restart the processing machine.
- ▷ Restart the PC.

9.2 Repairs

Repairs may be carried out only by trained, specialised personnel.

9.3 Warranty

The warranty complies with the statutory provisions. For more information, please refer to our General Terms of Business.

10 Environmental Protection

10.1 Packaging

In terms of packaging, the manufacturer participates in country-specific recycling systems, which ensure optimal recycling.

All packaging materials used are environmentally friendly and recyclable.

10.2 Waste equipment

Waste equipment contain materials that should be sorted for environmental-friendly recycling.

- ▷ After the service life has elapsed, dispose of the machine via the public disposal systems in an environmentally-friendly manner.

The plastic parts are labelled. Thus, they can be sorted for recycling or disposal.

10.3 Coolant/Lubricant (C/L)

- ▷ Dispose of the C/L in accordance with the applicable regulations and laws.

Proper disposal of the C/L must be documented.

11 Technical Data, Spare Parts and Accessories

Subject to changes.

Technical data of the processing machine

	Unit	Ceramill Matron	Ceramill Matron Core
Art. No.	–	181371	181370
Dimensions (D x W x H)	mm	650 x 940 x 850	650 x 574 x 850
Weight	kg	150	110
Electrical connections	V/Hz	100-230/50-60	
Laser power	W	220	
Fuse (time-delayed)	A	3.15 / 6.3	
Axes		5+ 3	5
Accuracy	µm	< 10	
Allowable temperature range (operation)	°C	18 - 30	
Max. operating height above NHN (standard elevation zero)	m	2000	
Max. relative humidity	%	80 (non-condensing)	
Spindle torque	Ncm	12	
Spindle speed	rpm	80000	
Collet diameter	mm	6 (optional 3)	
Sound pressure level	dB(A)	60	
Interfaces	–	USB / Ethernet / WLAN / Tank control / Extraction control	
Compressed-air connection	–	Dry, clean compressed air	
– Pressure	bar	6	
– Flow rate	l/min	50	
Air purity			
– Solid contaminants	–	Class 3; even better 5 µm for solids	
– Water content	–	Class 4; Max. pressure dew point +3 °C	
– Total oil content	–	Class 2; max. oil content 0.1 mg/m ³	

Tab. 2

	Unit	Ceramill Matron	Ceramill Matron Core
Extraction	–	Suitable for zirconium-oxide dusts; HEPA microfilter (captures 99.97 % at a particle size of 0.3 µm)	
_ Diameter of extraction connection	mm		38
_ Suction capacity	m ³ /h		240

Tab. 2

Technical data of tank

	Unit	Value
Art. No.	–	181375
Dimensions (D × W × H)	mm	500 × 200 × 400
Weight	kg	5.2
Capacity of C/L container	l	20
Mesh size of filter	μm	50
Interfaces	–	Pump unit

Tab. 3

Technical data of pump unit (CoolantFlow)

	Unit	Value
Art. No.	–	181376
Dimensions (D × W × H)	mm	440 × 150 × 370
Weight	kg	9.9
Max. flow rate of pump	l/min	3.8
Mesh size of filter	μm	5
Interfaces	–	Processing machine

Tab. 4

System components

Art. No.	Designation
181375	Ceramill Tank Matron
181376	Ceramill CoolantFlow (pump unit)
181377	Ceramill Cart Matron (cart)
179210	Service unit

Tab. 5

Accessory parts

Art. No.	Designation
181389	Auto-calibration kit Matron
181378	Ceramill SnapMag Matron
181251	Collet torque wrench Syc.3/6
179210	Service unit
181391	Collet cleaning set
180490	Collet, 3 mm
180539	Collet, 6 mm
180485	Chip tray Matron
180530	WLAN antenna
179256	Filter bag
181202	Filter cartridges
179257	Extraction hose with brush
178650-NTR	Coolant/lubricant concentrate
178652	pH-value test strips, 100 pce.
178653	Nitrite test, 100 pce. test strips
1156901	Fuse, processing machine (T 3.15 A / 250 V)
180310	Refractometer
179602	Torque wrench, 2 Nm, size 3 mm

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Amann Girrbach AG
Amann Girrbach GmbH