

airN GO

BEAUTY N HEALTH

The supra and perio double effect



POWDERS RANGE



Discover the

air **N** GO

THE FIRST RANGE OF POWDERS TO FOCUS ON THE PATIENT'S WELL-BEING

Ensure the Beauty and Health of your patients' teeth and the sustainability of their implants thanks to the new generation of prophylactic powders, specially designed to offer you increased gentleness and an outstanding comfort.

Increase possibilities for care and treatment with the new AIR-N-GO® formulas.

- AIR-N-GO “**CLASSIC**” and “**PEARL**” powders for **SUPRA**-gingival prophylaxis treatments.
- AIR-N-GO “**PERIO**” powder for **PERIO**dontal sub-gingival treatments.

Provide a moment of well-being to your patient!

50 powders

Natural Flavor
**TASTE
THE
DIFFERENCE**



They can

Dr Eric N., Bordeaux - France:

“With this new range of powders, I can envisage all types of treatment: prophylactic, aesthetic, paro-implantological! All these clinical applications can be envisaged with no special preparation. A possibility that I never had before! Whatever the type of patient who turns up at my surgery, I know that I've always got the powder that meets his needs ready to hand. ”

Stéphanie L., age 23, Marseille - France:

“I have a double problem: I'm hooked on caffeine and cigarettes! So it's not surprising: my smile was not very attractive, with stained, dull teeth. The results of this treatment are really impressive and visible! I no longer have the unsightly stains on my teeth that used to stop me from smiling normally. ”

Maya H., age 31, Brussels - Belgium:

“I love being able to choose the flavour I want! For once, I can have my say at the dentist's! I really liked the Raspberry - it was almost as if I was actually eating some! ”

Hugh M., age 34, Geneva - Switzerland:

“The flavours I was offered are incredible and amazing! I wanted to try them all! In fact, I got carried away and tried two: Peppermint and Raspberry, and I didn't have to wait long because the tanks of powders had already been prepared. I like Mint best, it's so refreshing! ”

Ludwig B., age 19, Stuttgart - Germany:

“Congratulations to my dentist! He's got what he needed to make treatment both efficient and pleasant! ”

not do without it

Carlos M., age 31, Barcelona - Spain:

“Like many people, I'm afraid of going to the dentist and I hate having my teeth polished! You have to admit that it hurts and isn't pleasant at all! But this time it was quite the opposite -completely different! During this session, I forgot my fear and felt no pain at all, though I'm very sensitive where my teeth are involved! It made me want to come back to my dentist! ”

Tina R., age 25, Berlin - Germany:

“The Cola flavour is extraordinary, almost a paradox at the dentist's surgery. With this little regression into infancy, I can almost forget I'm in your place! ”

Doris P., age 76, Zürich - Switzerland:

“At last, a Lemon flavour that tastes like lemon and doesn't leave a salty aftertaste! ”

Arnold HDL, age 64, Bordeaux - France:

“Doctor, your treatment is absolutely FANTASTIC!
It's as if I had new teeth! ”



"CLASSIC" AIR-N-GO

Maintenance and polishing: pleasant for all

This new concept of well-being relies on providing appropriate and tailored care to make daily prophylaxis easier.

The "CLASSIC" AIR-N-GO® powders offer a wide range of flavors, to discover depending on your patients' preferences.

The fine balance between flavor and sodium saccharin concentration removes any unpleasant salty taste and offers a feeling of freshness and immediate well.



Five fresh flavors based on 100 % natural aromas or essential oils:
Neutral, Cola, Raspberry, Peppermint, and Lemon.

100% NATURAL FRESHNESS!

AIR-N-GO powders

An unmatched comfort of the treatment

Thanks to the "CLASSIC" AIR-N-GO® powders, you will be able to show your patients, the benefits of prophylactic polishing.

The ultrafine grain size (76 µm) of the sodium bicarbonate and the less aggressive geometry of the crystals have been specifically studied to reduce the pain and bleeding caused by most competitors' powders.

The anti-clogging and controlled hydrophobic properties of these powders give them a dual action effect:


- in contact with water, micro particles shrink slightly: the polishing effect is much more efficient.
- the powder then dissolves gently, which avoids trauma to the delicate tissues and prevent risks of clogging of the AIR-N-GO device.

With the "CLASSIC" AIR-N-GO powders, you will immediately benefit from:

- a wide range of powders, made from natural components, with pleasant and surprising tastes which will delight your patients
- a thorough prophylaxis with instant and remarkable polishing results, even on difficult areas
- an efficient and gentle treatment that improves patient comfort: less bleeding and pain during treatment
- a steady spray stream for constant efficiency all day long
- a natural color recovery of teeth to strengthen the self-confidence of your patients
- a decreased risk of clogging




COMPETITOR

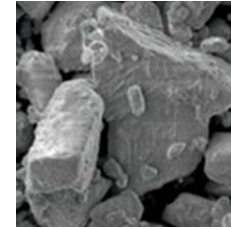

Competitor's powder:
layered particles with
aggressive geometry


x 1000

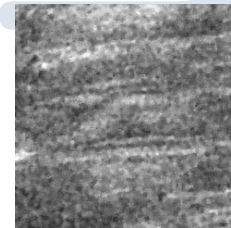

Enamel treated with
competitor's powder:
presence of numerous
scratches


x 300

air N GO

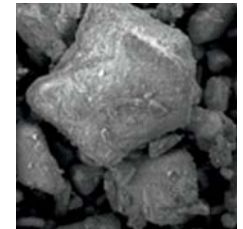

AIR-N-GO "CLASSIC" powder:
rounded geometry
particles

x 1000

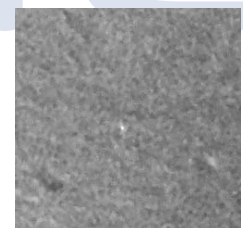

Enamel treated with
AIR-N-GO powder:
protected surface,
without impacts or scratches

x 300

Source : ICMCB_CNRS SEM photos

"PEARL" AIR-



Ultra-active polishing and protective powder

The less aggressive polishing effects of this powder, reveals the dazzling smile of your patients, without any pain.

The particular composition of the "PEARL" powder, in association with the AIR-N-GO® device, can restore the natural and pearly luster of the tooth enamel, in a smooth way.

Take advantage of a simple, ultra-fast, and remarkably effective treatment!

HIGH PERFORMANCE

AIR-N-GO powder


Efficiency and gentleness

Based on natural calcium carbonate, the “PEARL” powder, neutral tasting, consists of microspheres (55 µm) that are protective to delicate tissues. Its spherical structure allows easy contact with the tooth surface and the respect of the surrounding tissues.


Its total lack of sodium provides a naturally pleasant sensation in the mouth and is suitable for patients on a low-sodium diet.

Its unique properties give it many advantages:

- the enamel stays intact
- the lack of sodium makes the treatment applicable to the majority of patients
- the ultra-active polishing with respect for dental structures is painless for the patient
- the fast and secure removal of stains, plaque, extrinsic discoloration and stubborn deposits
- the tooth enamel gains a pearly appearance and a whiteness recovery
- the oral health is improved

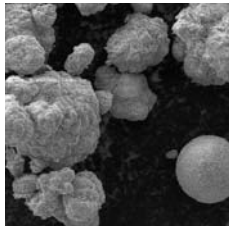

AIR-N-GO “PEARL” x 1000

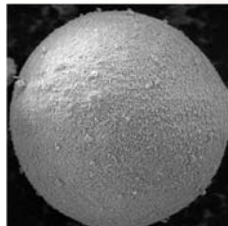

AIR-N-GO “PEARL” x 2500

"CLASSIC" and "PEARL"

No stress prophylaxis

AIR-N-GO® "CLASSIC" and "PEARL" powders have been specially developed for supra-gingival prophylaxis treatments. They are indicated for maintenance, caries prevention and improving oral hygiene of the patient.

The Beauty

In association with the AIR-N-GO device, in **SUPRA** function: these powders give cosmetic results thanks to a fast, effective and painless polishing that is gentle to the gingiva and teeth.

Protocol

POLISHING

To finish after scalling, to prepare for bleaching.

CLEANING

Interproximal areas, fissures and troughs, preparation of the thooth surface before etching. orthodontic brackets.

REMOVAL

Biofilm, plaque, stains and remnants ot temporary cementum.


clinical protocol

Power and preservation

The 120° and 90° AIR-N-GO® SUPRA nozzles associated with these powders ensure the ideal combination for any treatment of supra-gingival prophylaxis.

Protocol for using

- Direct the nozzle towards the tooth enamel allowing a distance of 3 to 5 mm.
- Vary the angle between the nozzle and tooth between 30 - 60 degrees.
- While administering treatment, make small circular movements over the tooth.
- Direct spray toward the occlusal edge.

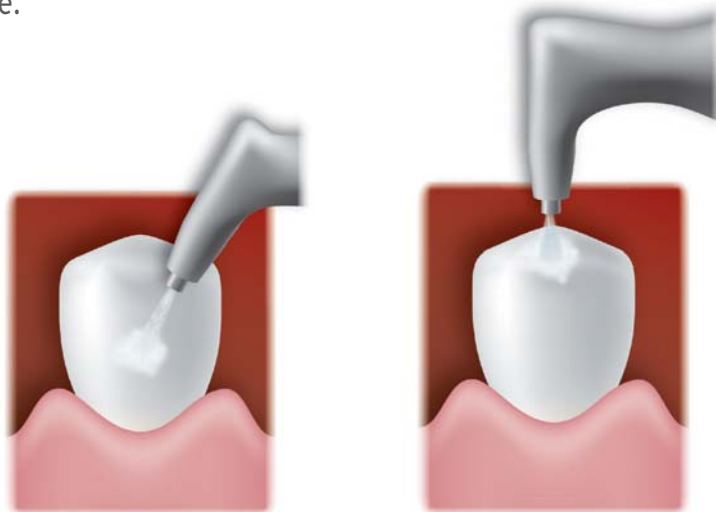

120° Standard

90° Optional

"PERIO" AIR-N-GO

Periodontal maintenance with exceptional results

The "PERIO" AIR-N-GO® powder, based on natural glycine, is used exclusively for sub-gingival treatment, in addition to initial therapy. Its efficacy has been proved by numerous studies and clinical trials*.

The low abrasiveness of the "PERIO" AIR-N-GO (25 µm) powder preserves the tooth enamel and the most delicate anatomical areas.

The Health

In association with the AIR-N-GO device, in the PERIO version: for a simple and efficient way to treat periodontal diseases and peri-implantitis sub-gingivally.



20 SECONDS per tooth or implant

* All references on our website: www.airngoconcept.com

O powder




Spectacular clinical actions:

- the prevention and treatment of periodontitis and peri-implantitis can be achieved in 20 seconds!**
- the sub-gingival plaque removal rapidly halt the spread of the biofilm.
- a sharp decrease of the inflammation and a reduction of even the deepest pocket are achieved in a few days, for periodontitis or peri-implantitis.
- the patient receives fast and painless treatment with guaranteed results in the treatment of periodontal pockets greater than 3 mm.

A technological revolution

The use of the PERIO nozzle combined with the “PERIO” powder allows the formation of a powerful and ultra-efficient spray on bacterial removal.

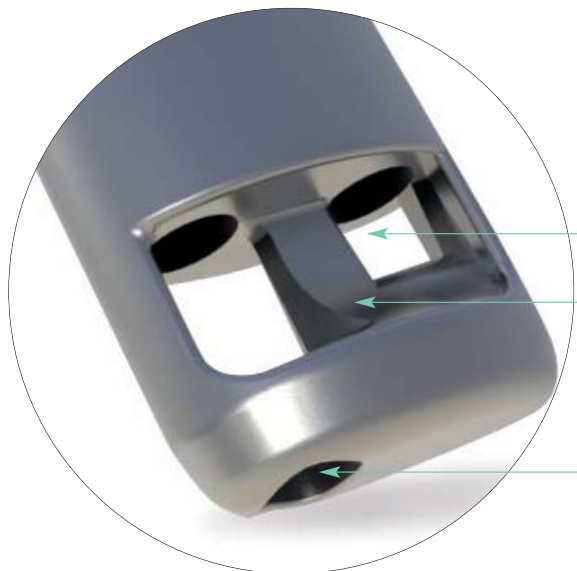

AN EXCLUSIVE SPRAY FORMATION PROCESS

Air, water and “PERIO” powder mixture chamber for a regular and constant spray.

Spray side deflector for a targeted cleaning action on the periodontal pockets.

“STOP POWDER” exit: for an air-water syringe action on demand.



"PERIO" clinical

A breakthrough for a revolutionary sub-gingival treatment

The "PERIO" AIR-N-GO® powder has been exclusively studied for the sub-gingival treatment of periodontal diseases and peri-implantitis. It can be used as a supplement in the initial therapy for pockets greater than 4 mm, when inflammation and / or bleeding have ceased.

Protocol

1ST STEP:

Probing, cleaning and debriding the pocket with SATELEC® ultrasonic tips/or manual instrumentation.



2ND STEP:


The AIR-N-GO "PERIO" powder is usually used seven to ten days after the 1st appointment, when inflammation and / or bleeding have ceased.

You must ensure that the process of attachment of collagen fibers on the root has not begun and that the tissues have minimal inflammation.

protocol

Protocol of using

- Use the surgical aspiration canula and maintain it near the treated area for cleaning the operatory area.
- Treat your patient as required.
- Each tooth or implant must be treated individually.
- Make sure there is no excessive bleeding during probing.
- Insert the nozzle to the bottom of the pocket exerting the same amount of pressure as with a periodontal probe.
- Activate the AIR-N-GO® “PERIO” device with your pedal and only when the nozzle is inside the pocket.
- Gently move the nozzle upwards moving it slowly from side to side without pressure for five seconds per side/surface.
- Leave AIR-N-GO activated even when you have gone past the top of the pocket. The gingival festoon must be treated slightly in order to obtain an optimal result.
- Repeat the operation on the surfaces of the tooth or implant to be treated.

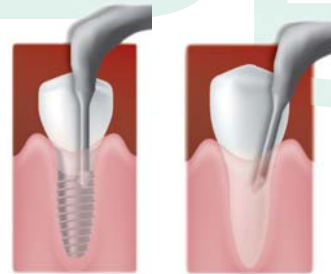

Treatment of peri-implantitis
Penetration of the PERIO nozzle
into a 5 mm pocket
Clinical case courtesy
of Dr E. Normand



Note: The “PERIO” AIR-N-GO powder can be used with the AIR-N-GO device in the SUPRA-gingival configuration for prevention and control periodontitis and perio-implantitis.

air N GO

BEAUTY N HEALTH

The supra and perio double effect



air N GO Powders

"CLASSIC"

"PEARL"

Box of 4 x 250 g transparent bottles.

Box of 50 x 20 g single use sticks.



air N GO Powder
PERIO

Box of 4 x 100 g
transparent bottles

Compatible devices



air MAX



prophy MAX



air N GO



air N GO
PERIO



air N GO

Find out more from: www.airngoconcept.com